Water Baths



Convincing quality features stand for a practiceproven range of 13 special Water Baths that guarantee reliable temperature examinations, independent of their field of application.





Quality built on Tradition

Users in research labs, standard and special labs for medicine, science and industry throughout the world have been profiting from the precision and dependability of our products for more than 45 years, all of which comply with valid European standards and bear the CE mark.

In the future, the success of a product spectrum of laboratory devices that meet the toughest demands on material, functioning and design will continue to be ensured by experience of every-day use in detail, ongoing technical advances and an excellent quality approach.

The phrase "Quality built on Tradition" encompasses more than just the sum of impressive product innovations. It is also an expression of our corporate policy, which includes a high degree of ready-and-waiting service as our primary customerfriendly service goal.

This applies to all of the GFL laboratory products that are produced exclusively at our plant in Germany - Deep Freezers, Shakers, Shaking Water Baths, Water Stills, Incubators or Water Baths alike.



Secure future

A vested quality demand in accordance with international standards is documented for all GFL laboratory products with the certification to DIN EN ISO 9001:2008, promoting more trust in the permanent high quality level of our products.

In addition to the continuous optimisation of ongoing production processes, the ISO obligation also calls for the commitment to quality awareness of our employees and the continuous development and rapid implementation of preventive measures to ensure quality assurance at a high level.

Valued as universal laboratory aids



In daily practice, GFL Water Baths have been valued for decades as a robust and universal aid for reliable temperature examinations in laboratories. They prove themselves daily in various fields of application – no matter whether they are used as

Incubation or Inactivation WaterBaths, as Steam Baths, as Water Baths for fume hoods, Multiple Water Baths or as a Tissue Float Bath. All GFL Water Baths are tested according to the German Equipment Safety Law and have the CE mark.

The accessories shown in this brochure are not included in the standard scope of supply of the units. The right of constructional modifications remains reserved. Special makes are available.



1002 - 1008

Successfully used

The range of Incubation and Inactivation Baths comprises seven models with five different sizes of 7 to 40 I capacity.

Water Bath 1005 is especially suitable for warming medical hot packs.

Water Baths 1012 and 1013 are equipped with a circulation system, thus guaranteeing an optimal temperature uniformity.



Specifications and features

- Short heating-up times.
- ▶ Temperature display and setting digitally via LED display, in 0.1 °C increments. Fast and exact setting, exact reproducibility of the set value.
- Electronic monitoring of the temperature regulator. In case of failure, the cause of the fault is shown on the display.
- Soft-touch keys with clear symbols. Two-finger operation prevents unintentional alterations of the set value.
- Interior (bath interior and heating element, cover frame, lid and perforated floor) made of stainless steel. Standard equipment: lid and perforated floor. Water Bath 1005 is supplied without perforated floor if the rack for hot packs 1923 (accessory) is ordered.

1005 Water Bath, 40 litres





Model Order No.	In Width mm	terior di Depth mm	mensions Working height • mm	Exter Width mm	ior dime Depth mm		Capacity approx. litres	Electrical connection*	Weight net kg	approx. gross kg	Packing volume approx m ³
1002	245	200	145	340	395	255	7	230 V / 5060 Hz / 1.0 kW	9,0	11	0.11
1003 1004	400 600	245 245	145 145	500 700	440 440	255 255	14 21	230 V / 5060 Hz / 1.5 kW 230 V / 5060 Hz / 1.5 kW	12,5 16.7	15 21	0.13 0.26
1005	410	296	315	510	490	445	40	230 V / 5060 Hz / 1.5 kW	20,7	27	0.30
1008	400	245	205	500	440	325	20	230 V / 5060 Hz / 1.5 kW	14,7	17	0.13
	• perfo	orated flo	or to upper rim					* Other v	oltages a	vailable o	n request


1012/1013

- Cover frame and lid have mirror finish. Lift-up, insulating lid with inner camber, no dripping back of condensate into the tubes.
- Corrosion-resistant outer housing, made of electrolytically galvanised sheet steel, powder-coated. Drain cock at the back.

Applications

All Water Baths models 1002 to 1013 are successfully used for incubations and inactivations of cultures, such as e.g. for warming bacteriological media, carrying out chemical reactions or thawing samples.

The circulation system of Water Baths 1012 and 1013 ensures optimum temperature uniformity throughout the whole bath. An electric motor with rotary magnet is flanged to the bath floor. Its torque is transmitted to a PTFE coated stirring magnet in the bath. The stirring magnet forces the water into the centre of the bath and then evenly back again.

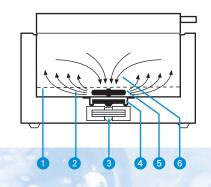

- 1 Perforated floor
- 2 Water circulation
- 3 Motor
- 4 Rotary magnet
- 5 Stirring magnet
- 6 Water circulation

Additional technical data of models 1002 - 1013

Temperature regulation: Microprocessor-controlled Temperature constancy: ± 0.1 °C temporal at 50 °C Temperature range: approx. 5 °C above ambient

temperature to +99.9 °C

with water level regulator 1919: approx. 3 $^{\circ}\text{C}$ above tap water

temperature to +99.9 °C

Over-temperature cut-out: electronic, 4 °C above set

temperature and electromechanical >130 °C



1013 Water Bath with circulation system, 14 litres

Model	In	terior d	imensions	Exter	ior dime	ensions	Capacity	Electrical connection*	Weight	approx.	Packing
Order	Width	Depth	Working height •	Width	Depth	Height	approx.		net	gross	volume
No.	mm	mm	mm	mm	mm	mm	litres		kg	kg	approx m ³
1012	245	200	145	340	395	325	7	230 V / 5060 Hz / 1.0 kW	10,6	12	0.11
1013	400	245	145	500	440	325	14	230 V / 5060 Hz / 1.5 kW	14,6	17	0.13
	◆ perforated floor to upper rim * Other voltages available on request										

Accessories



1002-1008



Test Tube Rack

made of stainless steel, with 20 openings of 18 mm \emptyset , for test tubes of 16/17 mm \emptyset , max. height 185 mm

Order No. 1920

Capacities:

 Water Bath 1002
 4 racks = 80 test tubes

 Water Bath 1003 / 1008
 8 racks = 160 test tubes

 Water Bath 1004
 12 racks = 240 test tubes

 Water Bath 1005
 9 racks = 180 test tubes



Test Tube Rack

made of stainless steel, with 5 openings of 31 mm Ø, max. height 185 mm

Order No. 1921

Capacities:

 Water Bath 1002
 4 racks = 20 test tubes

 Water Bath 1003 / 1008
 8 racks = 40 test tubes

 Water Bath 1004
 12 racks = 60 test tubes

 Water Bath 1005
 9 racks = 45 test tubes



Test Tube Rack

made of stainless steel, with 20 openings of 13 mm Ø, for test tubes of 12 mm Ø, max. height 185 mm

Order No. 1922

Capacities:

 Water Bath 1002
 4 racks = 80 test tubes

 Water Bath 1003 / 1008
 8 racks = 160 test tubes

 Water Bath 1004
 12 racks = 240 test tubes

 Water Bath 1005
 9 racks = 180 test tubes



Milk Bottle Rack

made of stainless steel, with 12 openings of 56 mm Ø, for baby milk bottles

Order No. 1942

Capacities:

Water Bath 1002 1 rack = 12 bottles
Water Bath 1003 / 1008 2 racks = 24 bottles
Water Bath 1004 3 racks = 36 bottles
Water Bath 1005 2 racks = 24 bottles




Adjustable Water Level Regulator

to keep the water level constant and to cool Water Baths 1002-1008 (approx. 3 °C above tap water temperature)

Order No. 1919



1005

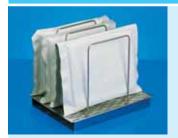

Rack for Water Bath 1005

made of stainless steel, especially suitable for warming medical hot packs.

Order No. 1923

Capacities:

Water Bath 1005 up to 8 hot packs

Accessories



1012/1013



Test Tube Rack

made of stainless steel, with 20 openings of 18 mm Ø, for test tubes of 16/17 mm Ø, max. height 185 mm

Order No. 1920

Capacities:

Water Bath 1012 4 racks = 80 test tubes Water Bath 1013 8 racks = 160 test tubes




Test Tube Rack

made of stainless steel, with 5 openings of 31 mm Ø, max. height 185 mm

Order No. 1921

Capacities:

Water Bath 1012 4 racks = 20 test tubes Water Bath 1013 8 racks = 40 test tubes

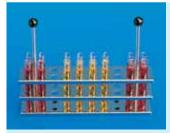

Test Tube Rack

made of stainless steel, with 20 openings of 13 mm Ø, for test tubes of 12 mm Ø, max. height 185 mm

Order No. 1922

Capacities:

Water Bath 1012 4 racks = 80 test tubes Water Bath 1013 8 racks = 160 test tubes



Milk Bottle Rack

made of stainless steel, with 12 openings of 56 mm Ø, for baby milk bottles

Order No. 1942

Capacities:

Water Bath 1012 1 rack = 12 bottles Water Bath 1013 2 racks = 24 bottles



Adjustable Water Level Regulator

to keep the water level constant and to cool Water Baths 1012-1013 (approx. 3 °C above tap water temperature)

Order No. 1919

1002-1013



Water Bath Protection Agent

ProAquaTop prevents the formation of algae, bacteria and mould. Biodegradable.

Please use biocides safely.
Always read labels and product information before use.

Highly efficient: only 1 ml per litre of water. A necessary change of water is indicated by fading of the blue colouration

Order No. 1910

Packing unit 1 bottle / Contents 200 ml.

Order No. 1911

Packing unit 3 bottles / Contents 3 x 200 ml

Order No. 1912

Packing unit 6 bottles / Contents 6 x 200 ml



Price available on request

please state number and diameter of openings

Lid with openings

made of stainless steel, with sets of rings made of heat-resistant plastic material. Number and diameter of openings can be chosen between 52 and 192 mm, in steps of 20 mm. Increases the variability of the Water Bath and reduces the loss of heat when placing large vessels into the Bath

Indispensable


For gentle steaming work with Erlenmeyer flasks, glasses, etc., GFL Water Baths models Steam Bath 1023 and Water Baths for fume hoods 1031 and 1032 have become indispensable in daily practice.

Accessories



Test Tube Rack made of stainless steel, with 100 openings of Ø 18 mm, for test tubes

Order No. 1933

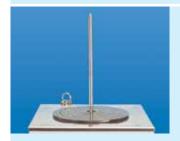

Support Rod made of stainless steel, length 316 mm, Ø 12 mm

GFL Water Baths convince through their varied ranges of application. For instance, the Steam Bath's square cover frame (W 265 mm x L 265 mm) is removable. Its 9-part set of rings, made of heat-resistant plastic material can be parted. Its diameter can be varied in approx. 20 mm steps (min. 32.5 / max. 173.5 mm).

Order No. 1985





Model Order No.			Exter Width mm	ior dime Depth mm	nsions Height mm	Capacity approx. litres	Electrical connection*	Weight net kg	approx. gross kg	Packing volume approx m ³	
1023	240 ♦ perfo	240 orated flo	120 or to upper rim	342	400	180	7	230 V / 5060 Hz / 1.0 kW * Other v	9 oltages a	12 vailable c	0.10 on request

1031/1032

Specifications and features

- ➤ Temperature range from approx. 5°C above ambient to boiling point, temperature regulation through tension thermostat. Heating element protected by over-temperature cut-out.
- ▶ 1023: bath, cover frame and perforated floor above the heating element are made of stainless steel.
- ▶ 1023: the adjustable water level regulator is situated at the back of the unit. It is part of the standard scope of supply, as well as the set of rings made of heat-resistant plastic material.
- ▶ 1023: outer housing made of electrolytically galvanised sheet steel, powder-coated.

- ▶ 1031/1032: bath, removable recessed lid with holes and perforated floor above the heating element are made of stainless steel.
- ▶ 1031/1032: the adjustable water level regulator at the right-hand side of the unit is part of the standard scope of supply.
- ▶ 1031/1032: the openings in the lid are covered by sets of rings made of heat-resistant plastic material.
- ➤ 1031/1032: four non-slip plastic feet ensure great stability.

Applications

Water Baths for fume hoods 1031 / 1032 are special baths that are preferably used for protected applications in fume hoods (digestories).



1031 Water Bath for fume hoods

Tec	hn	ical	ᅥᇰ	ta
166		IICa	ı ua	ιa

Model Order No.	Openings with sets of rings	Ø of the openings mm	Working height [♦] mm	Exter Width mm	ior dime Depth mm		Electrical connection*	Weight net kg	approx. gross kg	Packing volume approx. m ³
1031	6	91	100	430	300	155	230 V / 5060 Hz / 1.5 kW	6.1	7.4	0.05
1032	8 ♦ perfo	111 orated floor to	100 upper rim	670	300	155	230 V / 5060 Hz / 1.5 kW * Ot	8.2 her volta	10.2 ges availat	0.24 ble on request

1041/1042

Especially suitable

GFL Multiple Water Baths 1041 / 1042 are especially suitable for gentle steaming work with columns, Erlenmeyer flasks or beakers. Even after years of continuous operation they remain to be the most reliable laboratory aids.

Specifications and features


- Temperature range from approx. 5 °C above ambient to boiling point, temperature regulation through tension thermostat. Heating element protected by over-temperature cut-out.
- Bath interior, frame with openings and parted perforated floor above the heating element are made of stainless steel.
- The openings in the frame have a diameter of 130 mm and are covered by sets of rings made of heatresistant plastic material.
- ➤ A support rod of 600 mm length and a diameter of 12 mm, made of stainless steel, is situated behind each opening. It serves to fix the steaming vessels securely.

► The adjustable water level regulator is on the left-hand side of the bath. It is included in the standard scope of supply of the units.

➤ Outer housing made of electrolytically galvanised sheet steel, powdercoated.

Applications

Chemical, clinical and biological laboratories, both in research and industry, decide on our Multiple Water Baths whenever special quality and quantity requirements are to be fulfilled.

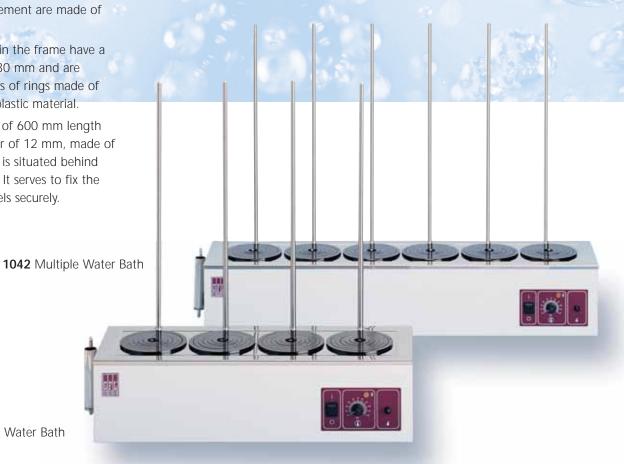

1041 Multiple Water Bath

Model Order No.	Openings with sets of rings	Ø of the openings mm	Working height [♦] mm	Exter Width mm	ior dime Depth mm	nsions Height mm	Electrical connection*	Weight net kg	approx. gross kg	Packing volume approx.
1041	4	131	90	682	232	190	230 V / 5060 Hz / 1.0 kW		15	0.12
1042	6 ◆ perfo	131 erated floor to	90 upper rim	982	232	190	230 V / 5060 Hz / 1.5 kW * Otl		20 ges availat	0.15 ble on request

Preferential choice

Special bath Tissue Float Bath model 1052 reliably fulfills everyday tasks in laboratories. Its exact temperature regulation ensures evenly stretched cuttings which are clearly visible in the black anodised bath. The cuttings are carefully dried on the warmed rim of the bath, which is also black anodised.

Specifications and features

- ► Exact temperature constancy of ±0.5° C.
- ➤ Temperature range from approx. 5 °C above ambient to approx. +80 °C.
- The low 100 mm height of the bath allows comfortable and safe working.
- Bath interior made of aluminium, black anodised. Housing made of aluminium, powder-coated, with knob and scale of the temperature regulator, mains switch and two pilot lamps (green/yellow, for mains and heating).
- ➤ The control thermometer with stainless steel holder is fixed on the rim of the bath.

Applications

Histological, chemical, clinical and bacteriological laboratories use this special bath to stretch and dry cut tissues.





Model Order No.	Temperatur range to	Interior dimensions Ø / H mm	Exterior dimensions Ø / H mm	Electrical connection*	Weigh net kg	t approx. gross kg	Packing volume approx m ³
1052	ca. +80 °C	200 / 60	280 / 100	230 V / 5060 Hz / 0.3 kW	2 * Other volta	5 iges availab	0.03 le on request



GFL Gesellschaft für Labortechnik mbH

P.O. Box 11 52 \cdot 30927 Burgwedel / Germany Schulze-Delitzsch-Strasse 4 \cdot 30938 Burgwedel / Germany Phone ++ 49 (0)5139 99 58 -0 \cdot Fax ++ 49 (0)5139 99 58 21

E-Mail: info@GFL.de · Internet: www.GFL.de